

GEM[®] PREMIER[™] 4000 with IQM[®]

Advanced, simple, and revolutionary.

Automation and innovation
at the touch of a button.



GEM[®] PREMIER[™] 4000

Now, automation is even better.

NEW

Our advanced PLUS technology enables:

- ▶ Faster touch-screen responsiveness
- ▶ Integrated wireless
- ▶ Remote service

PLUS 

Complements other key features:

iQM

Automates quality assurance

+

GEMweb Plus Custom Connectivity

Automates information management

+

Multi-use PAK

Automates maintenance

▼

Just press GO!

Automation and quality care go hand in hand.



Ultimate Flexibility

- ▶ Standardized platform for the lab and at the point-of-care
- ▶ Full analyte menu and integrated CO-Oximetry
- ▶ Accepts a syringe, capillary tube, or ampule
- ▶ Multiple menu and volume PAK configurations



MENU	PAK SIZE					Onboard Use-life
	75	150	300	450	600†	
Blood Gas, Hct, tHb, O ₂ Hb, HHb, COHb, MetHb, sO ₂ , Total Bili*	✓	✓	✓	✓	✓	30 days
Blood Gas, Electrolytes, Hct, tHb, O ₂ Hb, HHb, COHb, MetHb, sO ₂ , Total Bili*	✓	✓	✓	✓	✓	30 days
Blood Gas, Electrolytes, Glu, Lac, Hct, tHb, O ₂ Hb, HHb, COHb, MetHb, sO ₂ , Total Bili*	✓	✓	✓	✓	✓	30 days

Menus with: Creat, BUN, and measured tCO₂ are in development.

* Cartridges available with or without Total Bili.
 † Onboard use-life is 21 days.


Automated Quality Assurance

iQM: Intelligent Quality Management

The most accurate results, every time, for maximum efficiency and enhanced patient care.

Automatically monitors all components. Provides continuous error detection and correction.

iQM *Continuous Model*

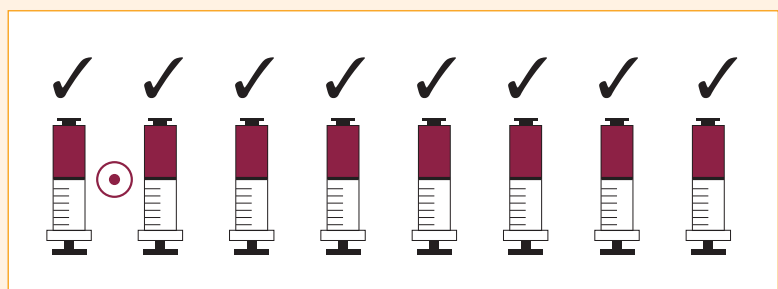


⊙ = Error automatically detected, corrected and documented immediately

✓ = Optimal test results ensured

0 —————> Continuous
Time

Traditional QC (Manual or Auto) *Eight-hour Model*

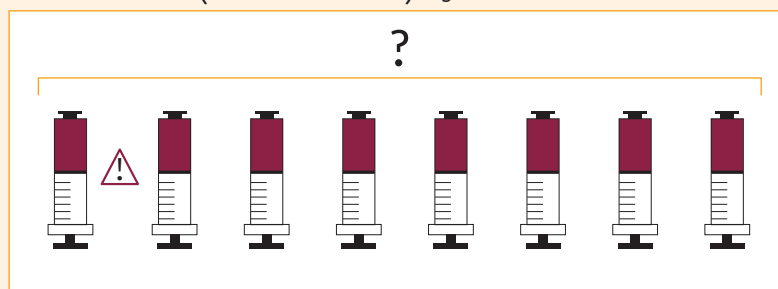

⚠ = Error undetected

? = Questionable test results (*all results from eight-hour period require review*)

0 —————> 8 hrs
QC IN Time QC OUT

Reduces Error Detection Time from Hours to Minutes^{1, 2}

	pH	pO ₂	pCO ₂	Na ⁺	K ⁺	Ca ⁺⁺	Glu	Lac	Hct
iQM	3 mins	3 mins	3 mins	17 mins	3 mins	3 mins	11 mins	6 mins	3 mins
Traditional QC (Manual or Auto)	≥8 hrs	≥8 hrs	≥8 hrs	≥8 hrs	≥8 hrs	≥8 hrs	≥8 hrs	≥8 hrs	≥8 hrs

Statistical presentation of an average error detection time with 95% confidence.

“iQM: a new standard for the future of QC.”


– James Westgard, PhD, developer of the ‘Westgard Rules’

References:

- Westgard JO, Fallon KD, Mansouri S. Validation of iQM Active Process Control Technology. *Point of Care, The Journal of Near-Patient Testing and Technology*, 2003: Vol. 2, No. 1.
- Toffaletti JG, McDonnell EH, Ramanathan LV, Tolnai J, Templin R, Pompa L. Validation of a quality assessment for a blood gas and electrolyte testing. *Clinica Chimica Acta*, 382 (2007) 65–70.

FDA-cleared in 2002 for GEM Premier 3000, and again in 2006 for GEM Premier 4000, **iQM is the only quality management system designed to provide continuous monitoring of the analytical process** with real-time, automatic error detection, correction and documentation of all corrective actions, replacing the use of traditional external quality controls.

Positive Impact on Staff Time

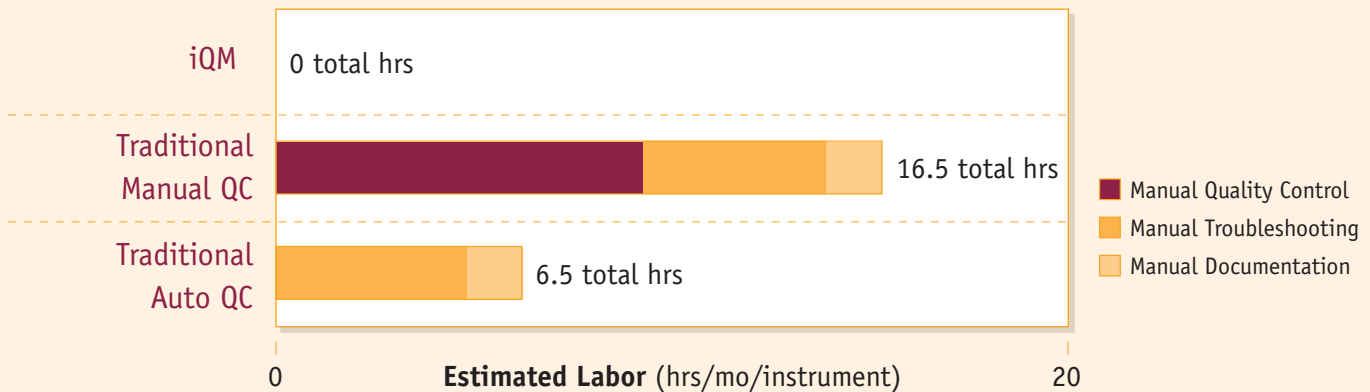

- ▶ Over 16 hours saved monthly
- ▶ Over 190 hours saved annually

This means:

In a hospital with 8 GEM Premier 4000 analyzers, iQM can save over 1,500 hours annually



IL's patented iQM and all analytical components are included inside a single-component, self-contained PAK.

Automated Information Management

GEMweb Plus Custom Connectivity

True web-based, bi-directional management from any networked analyzer or PC, customized for your hospital.

The ONLY connectivity software to:

- ▶ Provide a simple, single user-interface
- ▶ Allow total, system-wide control from any networked PC or GEM Premier 4000
- ▶ Automate and customize operator certification and management

AND:

- ▶ Allows automatic review of analyzer and PAK status
- ▶ Alerts operator to replace PAK
- ▶ Offers patient history with delta versus previous result with one keystroke
- ▶ Consolidates iQM reports, operator competency, audit logs and positive patient ID

Delivers automated quality assurance and regulatory compliance, while significantly reducing staff time.

Integrated wireless for seamless communication to the LIS or HIS.



GEMweb

Automated, web-based information management


- ▶ Single, graphic user-interface
- ▶ Bi-directional operator management
- ▶ Remote analyzer review and control
- ▶ Remote iQM review
- ▶ Demographic query



GEMweb Plus 100

Everything GEMweb offers AND:

- ▶ Remote global/area/analyzer system configuration
- ▶ Consolidated:
 - Analyzer status and control
 - Operator management
 - iQM reports
 - Patient information management
- ▶ System-wide sample lookup database
- ▶ HIS and LIS order receiving and processing

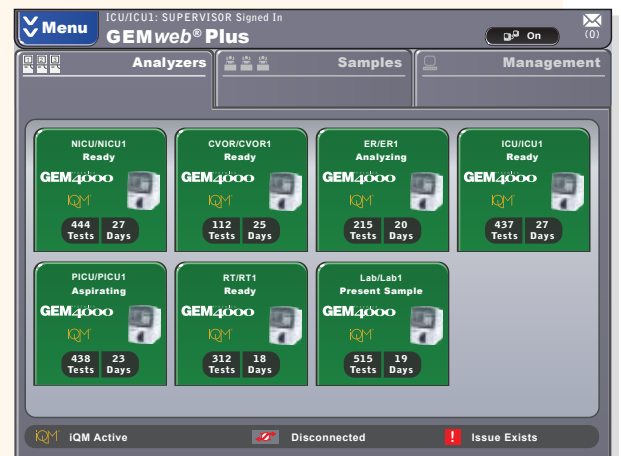

GEMweb Plus 200

Everything GEMweb Plus 100 offers AND:

- ▶ Operator competency ensured through automated and customizable certification, including exams
- ▶ Advanced regulatory compliance management
- ▶ Audit log creation



GEMweb Plus 100 & 200 Customizable Options

More functionality, more flexibility

- ▶ Receipt of ADT transmissions, including positive patient ID verification
- ▶ Automated and customizable point-of-care order creation (solicited and unsolicited)

Automated Maintenance

Saves significant staff time and reduces training and documentation requirements.

GEM PAK—The only single-component, multi-use cartridge on the market.

- ▶ Self-contained; includes all testing components: sensors, CO-Ox optical cell, lysing solution, Process Control Solutions, tubing, waste container and sampler
- ▶ Stored at room temperature
- ▶ Easy front-loading
- ▶ Replaced every 30 days

Automates the most labor- and skill-intensive processes.

- ▶ Zero maintenance—just replace the disposable PAK monthly; no additional cartridge handling required
- ▶ Only one PAK to inventory and manage, including all solutions, sensors, and quality control
- ▶ Immediately detects errors—no hands-on troubleshooting required
- ▶ Immediately corrects errors—no hands-on corrective actions required
- ▶ Regulatory compliance reports generated automatically



PAKs contain all components for patient testing and are ideal for high- and low-volume testing.

Automation at the Bedside

Bringing Simplicity to the Point-of-Care

Standardized lab-quality results and real-time connectivity for immediate decisions anywhere in the hospital.

CVOR

Immediate, precise testing

- ▶ CO-Oximetry, the most accurate method to measure Hemoglobin, is fully integrated, offering a complete test menu
- ▶ Measured Total Hemoglobin not affected by common operative variables, unlike conductive Hematocrit

ICU

Saves time and enhances patient care

- ▶ Provides complete oxygenation status
- ▶ Only cartridge-based system to provide Lactate testing in a single sample
- ▶ Comprehensive results in a single sample for fast clinical assessments and better outcomes
- ▶ Immediately recalls patient information, enhancing safety

NICU

Fast, reliable results for tiny patients

- ▶ The only cartridge-based system to offer Total Bilirubin
- ▶ Automatic Fetal Hemoglobin detection and correction
- ▶ Capillary-tube sampling (65 μ L)*
- ▶ Ideal for nitric oxide therapy

ED

Rapid, comprehensive testing†

- ▶ Only cartridge-based system to provide Lactate testing in a single sample—key to rapid implementation of sepsis protocol
- ▶ tHb for triaging trauma-related blood loss
- ▶ COHb, MetHb to assist in smoke-inhalation treatment
- ▶ Advanced GEMweb Plus connectivity—fully integrated, real-time information management



* BG, Electrolytes, Metabolites, Hct I † Basic Metabolic Panel in development: BUN/Creat and measured tCO₂

Simple and flexible, it's the automatic choice for your hospital.

The top choice of hospitals worldwide, the GEM Premier 4000 is advanced, simple, revolutionary and is fully automating critical care testing.

Automated Quality Assurance

iQM ensures test quality—24 hours a day, seven days a week, providing continuous error detection and correction.

Automated Information Management

GEMweb Plus Custom Connectivity offers true bi-directional communication anywhere—it automates information management and provides system-wide control from any networked PC or GEM Premier 4000.

Automated Maintenance

GEM PAKs contain all testing components, are replaced every 30 days and require no refrigeration. Combined with automated detection, correction and documentation, the most labor- and skill-intensive processes are eliminated, saving significant time.

Automation at the Bedside

Flexible cartridge test volumes and menu options, including a complete CO-Oximetry panel, and standardized, lab-quality results and real-time connectivity for immediate decisions at the POC.

Just press GO!



The GEM Family of Critical Care Testing Solutions

GEM Premier Analyzers	GEMweb Plus Custom Connectivity	GEM Complementary Products
<ul style="list-style-type: none">▶ GEM Premier 4000▶ GEM Premier 3500	<ul style="list-style-type: none">▶ GEMweb Plus 200▶ GEMweb Plus 100▶ GEMweb Plus Customizable Options▶ GEMweb	<ul style="list-style-type: none">▶ GEMdraw▶ GEM PCL Plus▶ GEM OPL

Every member of the GEM Family delivers fast, accurate results, in the lab or at the point-of-care, for the very best patient care.

Instrumentation Laboratory

Our Passion. Your Results.

Driven by our constant passion for innovation, IL has provided the most comprehensive and valuable diagnostic solutions for decades. From the launch of our first instrument, our integrated solutions continue to revolutionize critical care testing with more capabilities, more end-to-end automation and more flexibility.

Technical Specifications

Dimensions and Weight

Analyzer

H: 18 in, W: 12 in, D: 15 in, Wt: 44 lbs

PAK

H: 6.75 in, W: 10 in, D: 8 in, Wt: 8 lbs

Sample Volume

150 µL BG*/Hct/Lytes**/Glu/Lac/CO-Ox or any subset of the menu that includes CO-Ox

100 µL CO-Ox/Total Bili only

65 µL BG/Hct/Lytes/Glu/Lac

(micro mode) (capillary only)

*BG = pH, pCO₂, pO₂

**Lytes = Na⁺, K⁺, Ca⁺⁺, Cl⁻

Sample Type

Heparinized whole blood

Heparinized plasma (Total Bili)

Time to Results

All tests without CO-Ox: 70 seconds from sample introduction

All tests with CO-Ox: 95 seconds from sample introduction

Sample capacity: 75 – 600 tests

Measurement Methodology

Amperometric: pO₂, Glu, Lac, Creat[†]

Potentiometric: pH, pCO₂, Na⁺, K⁺, Ca⁺⁺, Cl⁻,
 BUN,[†] measured tCO₂[‡]

Conductivity: Hct

Optical measurement following chemical lysing of the whole blood sample:

CO-Ox, Total Bili

[†] In development. Not currently saleable.

Interface Protocols

ASTM or HL7 enables data transmission to a laboratory, hospital, or third-party information management system.

Measured Analytes

Analyte	Unit	Measured Range [‡]
pH	n/a	6.80 – 8.00
pCO ₂	mmHg	6 – 150
pO ₂	mmHg	5 – 800
Na ⁺	mmol/L	100 – 200
K ⁺	mmol/L	0.2 – 20.0
Ca ⁺⁺	mmol/L	0.10 – 5.00
Cl ⁻	mmol/L	40 – 170
Glu	mg/dL	4 – 750
Lac	mmol/L	0.3 – 20.0
Hct	%	15 – 75
tHb	g/dL	3.0 – 23.0
O ₂ Hb	%	-10.0 – 110.0
COHb	%	-10.0 – 110.0
MetHb	%	-10.0 – 110.0
HHb	%	-10.0 – 110.0
Total Bili	mg/dL	0.3 – 40
BUN [†]		
Creat [†]		
tCO ₂ [‡]		

[‡] The measured ranges reflect actual numeric value ranges the system can report.

Derived (calculated) Parameters

BE(B)	CvO ₂	paO ₂ /pAO ₂
BE(ecf)	p ₅₀	RI
tHb(c)	O ₂ cap	CcO ₂
Ca ⁺⁺ (7.4)	sO ₂	a-vDO ₂
Anion gap	sO ₂ (c)	Q _{sp} /Q _t (est)
P/F ratio	HCO ₃ ⁻ std	Q _{sp} /Q _t
pAO ₂	HCO ₃ ⁻ (c)	Hct(c)
CaO ₂	A-aDO ₂	

Werfen Group

Worldwide Locations

Corporate Headquarters

Barcelona, Spain
 Tel. +34-93-4010101
 www.werfengroup.com

Instrumentation Laboratory Headquarters

Bedford, MA
 Tel. +1-781-861-0710
 www.ilww.com

US, Canada, Latin America and South America

IL USA
 Bedford, MA
 Tel. +1-781-861-0710
 www.ilus.com

Werfen Medical Ltd. Brazil
 São Paulo
 Tel. +55-11-4154-3337

IL Canada
 Richmond Hill, ON
 Tel. +1-800-552-2025
 x6115

IL Mexico
 Col. Granada
 Tel. +52-55-5262-1760
 www.il-mexico.com.mx

Izasa Uruguay
 Montevideo
 Tel. +59-82-4818133

Pacific

Werfen Medical China
 Shanghai
 Tel. +86-21-32100745

Beijing
 Tel. +86-10-59756032

Werfen Hong Kong
 Hong Kong
 Tel. +852-27927773

IL India
 New Delhi
 Tel. +91-11-25191251

IL Japan
 Minato-ku, Tokyo
 Tel. +81-3-5419-1301

Werfen Medical Korea
 Seoul
 Tel. +82-2-571-9246
 www.werfenmedical.com

Europe, Middle East, Africa

IL Austria
 Vienna
 Tel. +43-1-2565800-0

IL Belgium
 Zaventem
 Tel. +32-2-7252052
 www.il-be.com

Comesa Czech
 Prague
 Tel. +420-2-7816047

IL France
 Le Pré Saint Gervais
 Tel. +33(0)1-82-308600
 www.il-france.fr

IL Germany
 Munich
 Tel. +49-89-909070
 www.il-ger.de

IL The Netherlands
 Breda
 Tel. +31(0)-76-5480100
 www.il-nl.com

Comesa Hungary
 Budapest
 Tel. +36-1-4392910 or 11

IL Italy
 Milan
 Tel. +39-02-25221
 www.il-italia.it

IL Lithuania
 Kaunas
 Tel. +370-37-313157

Comesa Poland
 Warsaw
 Tel. +48-22-3361800

Izasa Portugal
 Carnaxide
 Tel. +351-21-4247300
 www.izasa.com

IL Russia
 Moscow
 Tel. +7-495-1244559

Izasa Spain
 Barcelona
 Tel. +34-93-4010101
 www.izasa.com

IL UK
 Warrington, Cheshire
 Tel. +44-1925-81-0141
 www.il-uk.com

Instrumentation Laboratory, Comesa and Izasa are companies of Werfen Group.

For more information, call **1.800.955.9525**
 or visit **www.ilus.com**.

The Instrumentation Laboratory logo, GEM, Premier, iQM, GEMweb, and GEM PCL are trademarks of Instrumentation Laboratory Company and/or one of its subsidiaries or parent companies, and may be registered in the United States Patent and Trademark Office and in other jurisdictions. All other product names, company names, marks, logos, and symbols are trademarks of their respective owners.

©2012 Instrumentation Laboratory. All rights reserved.

